



ARTICLE


Rates of Recovery of *Enterobacter sakazakii* (*Cronobacter* spp.) from Powdered Infant Formula Using Both a Chromogenic Agar and Real-Time PCR : A Preliminary Study

Kwang-Young Song^{1,2†*}, Kun-Ho Seo^{1†}, and Jung-Whan Chon³

¹Center for One Health and College of Veterinary Medicine, Konkuk University, Seoul, Korea

²Department of Zoo Engineering, International University of Ulaanbaatar (IUU), Ulaanbaatar, Mongolia

³Department of Pet Total Care, Kyung-in Women's University, Incheon, Korea

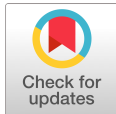

Received: September 3, 2021
Revised: September 19, 2021
Accepted: September 24, 2021

†These authors contributed equally to this study.

*Corresponding author :
Kwang-Young Song
Center for One Health and College of Veterinary Medicine, Konkuk University, Seoul, Korea
Tel : +82-2-450-4121
Fax : +82-2-3436-4128
E-mail : drkysong@gmail.com

Copyright © 2021 Korean Society of Dairy Science and Biotechnology.
This is an Open Access article distributed under the terms of the Creative Commons Attribution Non-Commercial License (<http://creativecommons.org/licenses/by-nc/3.0>) which permits unrestricted non-commercial use, distribution, and reproduction in any medium, provided the original work is properly cited.

ORCID

Kwang-Young Song
<https://orcid.org/0000-0002-5619-8381>
Kun-Ho Seo
<https://orcid.org/0000-0001-5720-0538>
Jung-Whan Chon
<https://orcid.org/0000-0003-0758-6115>

Abstract

Although the number of incidences of illness caused by ingestion of the bacterial pathogen *Enterobacter sakazakii* (*Cronobacter* spp.) has dramatically declined, there remains a need for a robust isolation method to recover this microbe from powdered infant formula (PIF). The current method described in the FDA's Bacteriological Analytical Manual requires multiple steps, and 3-4+ days for complete analysis of PIF isolated *E. sakazakii* (*Cronobacter* spp.). We describe a bacteriological method including a one-step enrichment followed by plating on chromogenic agar for presumptive identification of *E. sakazakii* (*Cronobacter* spp.). Suspected colonies are confirmed by either biochemical analyses, or a Real-Time PCR-based assay. Using this method, *E. sakazakii* (*Cronobacter* spp.) in PIF can be isolated and identified within one day (24 hours).

Keywords

Enterobacter sakazakii (*Cronobacter* spp.), enrichment, chromogenic agar, real-time PCR, Rapid ID 32E

Introduction

On April 12, 2002, the FDA had alerted health care professional about the risk of *Enterobacter sakazakii* (*Cronobacter* spp.) infection in hospitalized newborn infants, particularly premature infants or other immuno-compromised infants who fed powdered infant formulas [1-5]. This alert issued by FDA was targeted to concerns for immuno-compromised infants in hospital settings [1,3-5]. The conventional methods for isolating and identifying *E. sakazakii* (*Cronobacter* spp.) from infant formula required enrichment culture for about 48 h, and then subculture to selective agar followed by phenotypic identification [1-3]. Namely, these procedures took up at least 5 days to obtain a result [1-3]. Therefore, in order to isolate and identify *E. sakazakii* (*Cronobacter* spp.) rapidly, various studies have been conducted, and various attempts to inhibit the growth of *E. sakazakii* (*Cronobacter* spp.) had also been conducted [1,5-7]. Fortunately, Real-Time PCR could provide the rapid-accurate-quantitative analysis for confirming various food-borne pathogens [8-10]. The format was developed for detecting *E. sakazakii* (*Cronobacter* spp.) in dried infant formula using by applying the fluorogenic 5' nuclease assay (TaqMan[®]) with ABI Prism 7000 [10]. Seo and Brackett [10] firstly made the format of probe and primers based on the partial macromolecular synthesis operon

of *E. sakazakii* (*Cronobacter* spp.) such as the rpsU gene 3' end and the primase (dnaG) gene 5' end [10]. Also, in the Real-Time PCR instrument, the fluorescence measurement is performed in real time and automatically analyzed by the SDS software installed in the sequence detector [8-10]. Therefore, using a Real-Time PCR instrument, results could be quickly confirmed within 2 hours [8-10].

Also, various chromogenic agars have been developed to detect *E. sakazakii* (*Cronobacter* spp.). To date, various companies have developed chromogenic agar products for detection of *E. sakazakii* (*Cronobacter* spp.), but it is urgently needed to develop a method to solve the problem of false positives and false negatives [11-13]. Among the various chromogenic agars, Druggan-Forsythe-Iversen (DFI, UK) was first developed for the purpose of detecting *E. sakazakii* (*Cronobacter* spp.), and is still widely used.

Hence, this study developed and compared Real-Time PCR assays based on the macromolecular synthesis (MMS) gene, for the detection of *E. sakazakii* (*Cronobacter* spp.) using the ABI Prism 7000 SDS (Applied Biosystems, USA) instrument with a selective chromogenic agar within 24 hours. And also real-time PCR approaches were compared with the conventional culture method for the detection of *E. sakazakii* (*Cronobacter* spp.) in dried infant formula.

Materials and Methods

1. Bacterial strains

All strains (*E. sakazakii* [*Cronobacter* spp.] and non-*E. sakazakii* [non-*Cronobacter* spp.]) were provided by Division of Microbiology, CFSAN, FDA, and these bacterial strains were cultured in tryptic soy broth (TSB; Becton Dickinson, USA) at 37°C. To obtain an initial low cell number, 10-fold serial dilutions were performed in phosphate-buffered saline (PBS, pH 7.2). Viable cell numbers of these strains were determined by direct plating to DFI (Chromogenic *E. sakazakii* [*Cronobacter* spp.] agar, CM1055, Oxoid) following incubation at 37°C for 24 hours. Also, these strains were examined for its response to stress condition and recovery on selective media as DFI.

2. Spiking of the dried infant formula and limit of detection

Various types of dried infant formula were purchased from a grocery store in College Park, USA. Dried infant formula were artificially inoculated at several different target levels, and the inoculated test portions were submitted to the one-enrichment method at Fig. 1, using by DFI as selective media and real-time PCR. Assays were repeated 3 times for all contamination levels including for the blank (negative control). Detection limit was assessed for *E. sakazakii* (*Cronobacter* spp.) N6 strain. This stain was selected to present the typical yellow pigmentation on TSA, the green color on DFI, and positive on Real-Time PCR.

3. Confirmation by real-time PCR and rapid ID 32 E

A real-time PCR assay developed by Seo and Brackett [10] targeting the dnaG gene

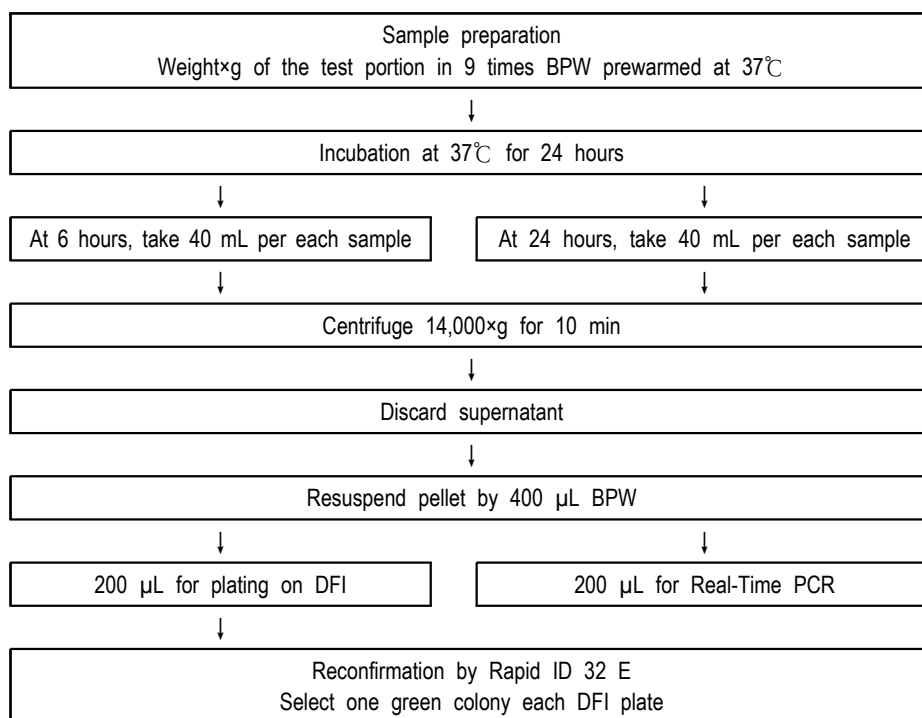


Fig. 1. The flow scheme of one-time enrichment method for detecting *E. sakazakii* (*Cronobacter* spp.) using by selective media as DFI and real-time PCR in dried infant formula. DFI, Druggan-Forsythe-Iversen.

on the MMS operon of *E. sakazakii* (*Cronobacter* spp.) was modified and used to confirm identification of each sample. Two sets of each sample (200 µL) from one-time enrichment broth were centrifuged at 14,000×g for 10 min.

Negative control (filter sterilized water) and positive control (1 mL of *E. sakazakii* (*Cronobacter* spp.) with over 10^8 CFU/mL) were included in every assay. The pellets of cell were resuspended in PreMan Ultra (Applied Biosystems) and then placed in the boiling water (about 100°C) during ten min. Hence, samples were cooled for 2 min at ambient temperature and also centrifuged at 14,000×g for 10 min again.

Forty-five microliter of PCR mix was combined by 5 µL of forward primer, 5 µL of reverse primer, 5 µL of TaqMan probe, 5 µL of filter sterilized water, and 25 µL TaqMan Universal PCR Master Mix in a 96-well microwell plate, respectively. And then, 5 µL of the supernatant fluids was transferred to the prepared PCR mix. Next, the 96-well microwell plates were sealed with optical adhesive covers and were analyzed by ABI Prism 7000 SDS (Applied Biosystems). The reaction was run at 50°C for 2 min and then 95°C for 10 min, followed by 40 cycles of 95°C for 15 s and 60°C for 60 s.

Furthermore, it was also investigated whether the ABI Prism 7000 SDS format designed for detecting *E. sakazakii* (*Cronobacter* spp.) was applicable to real-time PCR devices such as rotor gene, smart cycler, and light cycler produced by other companies.

Also, the presumptive colony of *E. sakazakii* (*Cronobacter* spp.) shown green color was taken from each DFI agar plate and then were biochemically identified by the Rapid

ID 32 E strips (bioMérieux), which were incubated at 37°C for approximately 4 hours.

4. Competition experiments with other Enterbacteriaceae and *Salmonella* spp.

Varied cell numbers of *E. sakazakii* (*Cronobacter* spp.) along with defined numbers of either *Salmonella* Enteritidis or *Enterobacter cloacae* were spiked simultaneously into (ex. 25)g of dried infant formula in 9 times (ex. 25×9 times) mL buffered peptone water (BPW), and then followed by one-time enrichment method (Fig. 1).

5. Statistical analysis

Two replicate experiments were performed for each bacterium studied. Data were analyzed using the statistical program, and Duncan’s multiple range test was used to separate means using a level of significance of $p < 0.05$.

Results & Discussion

The purpose of this study was to detect *E. sakazakii* (*Cronobacter* spp.) within 24 hours. The first is the result of a comparative experiment on the most suitable enrichment broth for *E. sakazakii* (*Cronobacter* spp.). First, mLST, universal pre-enrichment broth, water including 0.85% NaCl, buffered peptone water, widely used as an enrichment medium for detecting microorganisms, and R&F enrichment broth sold as a selective enrichment medium for *E. sakazakii* (*Cronobacter* spp.) were compared respectively (Table 1).

Growth of *E. sakazakii* (*Cronobacter* spp.) was confirmed through DFI, an *E. sakazakii* (*Cronobacter* spp.) selective chromogenic agar, while culturing for 0 hours, 6 hours, and 18 hours. mLST, R&F enrichment broth, universal pre-enrichment broth, water including 0.85% NaCl, and buffered peptone water were all detected at 18 hours, and only R&F enrichment broth, universal pre-enrichment broth, water including 0.85% NaCl, and buffered peptone water were detected at 6 hours. As a result, buffered peptone water showed the best results (Table 1).

In addition, an experiment was conducted in which vancomycin, an antibiotic capable of inhibiting Gram-positive bacteria, was added to various enrichment broth. As a result, there was no statistically significant difference between the absence of

Table 1. Comparison of detection when inoculated with less than 1.0 cfu/mL of *E. sakazakii* (*Cronobacter* spp.) N6 in various enrichment both

Incubation time	Type of enrichment media					Colony color on DFI
	M	R	U	W	B	
0 hour	-	-	-	-	-	-
6 hours	-	2.65	2.14	3.1	3.4	Green ¹⁾
18 hours	6.7	9.04	9.41	9.47	9.59	Green

(-) is no grow.

¹⁾ Green colony color recognizes as *E. sakazakii* (*Cronobacter* spp.) on DFI.

M, mLST; R, R&F enrichment broth; U, universal pre-enrichment broth; W, water including 0.85% NaCl; B, buffered peptone water; DFI, Druggan-Forsythe-Iversen.

vancomycin and the addition of vancomycin (final concentration is 8 µg/mL) (data not shown).

Next, the inclusivity and exclusivity were compared using various methods used to identify *E. sakazakii* (*Cronobacter* spp.) (Table 2). Summarizing the results obtained in this study, nutrient agar did not show any differentiation because all microorganisms were able to grow. However, since DFI is a chromogenic agar that can detect *E. sakazakii* (*Cronobacter* spp.), it was able to specifically observe green colony. Also, Ripid ID 32E was able to distinguish *E. sakazakii* (*Cronobacter* spp.). Furthermore, real-time PCR was able to confirm the presence of *E. sakazakii* (*Cronobacter* spp.) in a short time (Table 2).

Therefore, to summarize the results of this experiment, the optimal combination for detecting *E. sakazakii* (*Cronobacter* spp.) is to check whether positive or negative by the color of the colony with DFI agar, and to check whether it is *E. sakazakii* (*Cronobacter* spp.) or not by using real-time PCR at the same time. Then, if necessary, additional confirmation of green colonies on DFI agar by Rapid ID 32E method is considered to be able to accurately identify *E. sakazakii* (*Cronobacter* spp.).

In order to detect *E. sakazakii* (*Cronobacter* spp.) within 24 hours, the method of determining whether positive or negative for *E. sakazakii* (*Cronobacter* spp.) is positive or negative with various newly developed devices such as real-time PCR in the enrichment medium for 6 hours after starting the enrichment culture broth is considered very important. In particular, securing a large number of bacteria is the biggest challenge in order to confirm positive or negative of *E. sakazakii* (*Cronobacter* spp.) in the enriched culture broth for 6 hours. Therefore, the centrifugation technique, which is a method commonly used in *Escherichia coli* testing, was applied to the detection of *E. sakazakii* (*Cronobacter* spp.) [14]. Therefore, the enrichment culture was started for the samples to which various concentrations of *E. sakazakii* (*Cronobacter* spp.) were added, and the enriched culture broth was centrifuged at 6 hours and 24 hours, respectively, and those without centrifugation were compared with each other (Table 3).

Table 2. Inclusivity and exclusivity as confirmation of *E. sakazakii* (*Cronobacter* spp.)

Strain tested in study	Confirmation method				
	DFI ¹⁾	Nutrient agar	Rapid ID 32E	Real-Time PCR	
<i>Pantoea agglomerase</i>	Yellow	Yellow	Non- <i>E.sakazakii</i>	Negative	
<i>Salmonella</i>	Enteritidis 18088	Black	Yellow	Non- <i>E.sakazakii</i>	Negative
	Enteritidis 18673	Black	Yellow	Non- <i>E.sakazakii</i>	Negative
	Typhimurium LT2	Black	Yellow	Non- <i>E.sakazakii</i>	Negative
<i>Enterobacter</i>	<i>gergorivae</i>	Yellow	Yellow	Non- <i>E.sakazakii</i>	Negative
	<i>dissolvens</i> E10	Yellow	Yellow	Non- <i>E.sakazakii</i>	Negative
	<i>helveticus</i> E912	Green	Yellow	Non- <i>E.sakazakii</i>	Negative
	<i>helveticus</i> E440	Weak green	Yellow	Non- <i>E.sakazakii</i>	Negative
	<i>vulneris</i>	Green	Yellow	Non- <i>E.sakazakii</i>	Negative
	<i>sakazakii</i> N6	Green	Yellow	<i>E. sakazakii</i>	True positive
<i>sakazakii</i> 626	Weak green	Yellow	<i>E. sakazakii</i>	True positive	

¹⁾ Druggan-Forsythe-Iversen (DFI, Oxoid) as chromogenic agar was first developed for the purpose of detecting *E. sakazakii* (*Cronobacter* spp.) which was shown green colony color on DFI.

Table 3. Effect of detecting *E. sakazakii* (*Cronobacter* spp.) N6 between spin down and non-spin down samples (unit: log cfu/mL)

Inoculation concentration	<i>E. sakazakii</i> (<i>Cronobacter</i> spp.) N6			
	Spin down ¹⁾		Non-spin down as control	
	6 h	24 h	6 h	24 h
10 cfu/25 g	3.63–3.76	Over 9	2.30–3.00	Over 9
1 cfu/25 g	2.79–3.02	Over 9	1.60–1.90	Over 9
0 cfu/25 g	-	-	-	-

¹⁾ The volume of spin down is 10 mL.

The method of centrifuging the enriched culture broth for both 6 and 24 hours showed a good result, which was very helpful in determining *E. sakazakii* (*Cronobacter* spp.), as the number of bacteria was detected with a statistically significant difference ($p < 0.05$).

Next, in order to secure more bacterial counts of *E. sakazakii* (*Cronobacter* spp.), we compared the results of bacterial counts in DIF, a selective chromogenic agar, and real-time PCR for *E. sakazakii* (*Cronobacter* spp.), when the centrifuged amount of the enriched culture medium was expanded to 40 mL (Table 4). As a result, as the amount of centrifugation of the enrichment culture increased, better results were obtained (Table 4). In addition, according to previous study, the real-time PCR results were compared when the centrifuged amounts of the enriched culture medium were 1 mL, 10 mL, and 40 mL at 6 hours and 24 hours, respectively. As expected, the best results were obtained when the amount of centrifugation was 40 mL (Data not shown).

The real-time PCR used in this study was manufactured based on the ABI Prism 7000 SDS. Therefore, we tested whether this format of real-time PCR, which was made to identify *E. sakazakii* (*Cronobacter* spp.), can be applied to other real-time PCR devices (Table 5). This is because each laboratory has real-time PCR devices produced by various companies.

A total of four real-time PCR devices compared in this study are ABI Prism 7000 SDS, rotor gene, smart cycler, and light cycler. As shown in Table 5, the format of ABI Prism 7000 SDS, which can discriminate positive and negative for *E. sakazakii* (*Cronobacter* spp.), was directly applied to rotor gene and smart cycler to accurately discriminate positive and negative. However, it was shown that it was not applied in the case of light

Table 4. Detection limit of *E. sakazakii* (*Cronobacter* spp.) N6 between DFI and real-time PCR after 40 mL sample spin down

Inoculation concentration	<i>E. sakazakii</i> (<i>Cronobacter</i> spp.) N6			
	40 mL spin down (log cfu/mL)		P/T ¹⁾	
	6 h	24 h	6 h	24 h
56 cfu/30 g	5.78	9.83	2/2	2/2
7 cfu/30 g	5.22	10	2/2	2/2
2 cfu/30 g	5.05	10	2/2	2/2
0 cfu/30 g	-	-	0/2	0/2

¹⁾ Positive number per total number analyzed by real-time PCR. DFI, Druggan-Forsythe-Iversen.

Table 5. Comparison of the Ct value of *E. sakazakii* (*Cronobacter* spp.) and non-*E. sakazakii* (non-*Cronobacter* spp.) using by various real-time PCR machine

Type of real-time PCR machine	List of stains tested in this study (Ct value)										Volume of format
	<i>Escherichia vulneris</i>	<i>Enterobacter cloacae</i>	<i>E. sakazakii</i> (<i>Cronobacter</i> spp.)				<i>Salmonella</i> Enteritidis	Water	Positive control	Negative control	
			ES626	N6	FSM282	N6					
ABI Prism 7000 SDS	Undet ¹⁾	Undet	24.88	25.15	20.90	24.05	Undet	Undet	22.49	Undet	50 µL
Rotor-gene	Undet	Undet	17.03	23.59	15.90	17.25	Undet	Undet	15.76	Undet	25 µL
Smart cyclor	Undet	Undet	18.33	23.23	17.22	19.55	Undet	Undet	17.59	Undet	25 µL
Light cyclor	Undet	Undet	Undet	Undet	Undet	Undet	Undet	Undet	Undet	Undet	20 µL

¹⁾ undetection.

cyclor (Table 5).

In conclusion, this study established a method using a selective chromogenic agar and real-time PCR simultaneously to detect *E. sakazakii* (*Cronobacter* spp.) within 24 hours. An enrichment culture broth capable of selectively enriching *E. sakazakii* (*Cronobacter* spp.) was investigated. In addition, by applying the centrifugation method, *E. sakazakii* (*Cronobacter* spp.) was detected even if the incubation time was 6 hours. In order to further improve this method, it is necessary to set the number of selective chromogenic agar to at least two. Additional rapid & accurate detection methods should also be developed. In addition, further studies are needed to increase the sensitivity of real-time PCR for confirming *E. sakazakii* (*Cronobacter* spp.). Furthermore, it is necessary to actively use upgraded biochemical methods (for example, Matrix-Assisted Laser Desorption Ionization-Time-of-Flight Mass Spectrometer, etc) that are currently being developed in various ways [15]. Further research is needed so that the method developed in this study can be used in the detection of other food-borne poisoning bacteria.

Conflict of Interest

The authors declare no potential conflict of interest.

Acknowledgements

This work was supported by Korea Institute of Planning and Evaluation for Technology in Food, Agriculture, Forestry (IPET) through Useful Agricultural Life Resources Industry Technology Development Program, funded by Ministry of Agriculture, Food and Rural Affairs (MAFRA)(121051021HD020) and also by the Konkuk University Researcher Fund, 2021.

References

1. Pina-Pérez MC, Rodrigo D, Martínez A. Nonthermal inactivation of *Cronobacter sakazakii* in infant formula milk: a review. *Crit Rev Food Sci Nutr.* 2016;56:1620-1629.

2. Proudly I. *Cronobacter* spp. (“*Enterobacter sakazakii*” sensu lato): implication dans la contamination des préparations en poudre pour nourrissons et enfants en bas âge. *Can J Microbiol.* 2009;55:473-500.
3. Chon JW, Seo KH, Yim JH, Bae D, Kim B, Kim TJ, et al. Evaluation of selective media containing iron source and alpha-glucosidase substrates for *enterobacter sakazakii* (*Cronobacter* spp.) detection. *J Dairy Sci Biotechnol.* 2021;39:9-19.
4. Nazarowec-White M, Farber JM. *Enterobacter sakazakii*: a review. *Int J Food Microbiol.* 1997;34:103-113.
5. Chauhan R, Singh N, Pal GK, Goel G. Trending biocontrol strategies against *Cronobacter sakazakii*: a recent updated review. *Food Res Int.* 2020;137:109385.
6. Friedemann M. *Enterobacter sakazakii* in food and beverages (other than infant formula and milk powder). *Int J Food Microbiol.* 2007;116:1-10.
7. Henry M, Fouladkhah A. Outbreak history, biofilm formation, and preventive measures for control of *Cronobacter sakazakii* in infant formula and infant care settings. *Microorganisms.* 2019;7:77.
8. Zemanick ET, Wagner BD, Sagel SD, Stevens MJ, Accurso FJ, Harris JK. Reliability of quantitative real-time PCR for bacterial detection in cystic fibrosis airway specimens. *PLOS One.* 2010;5:e15101.
9. Catarame TMG, O’Hanlon KA, McDowell DA, Blair IS, Duffy G. Comparison of a real-time polymerase chain reaction assay with a culture method for the detection of *Salmonella* in retail meat samples. *J Food Saf.* 2006;26:1-15.
10. Seo KH, Brackett RE. Rapid, specific detection of *Enterobacter sakazakii* in infant formula using a real-time PCR assay. *J Food Prot.* 2005;68:59-63.
11. Gurtler JB, Beuchat LR. Performance of media for recovering stressed cells of *Enterobacter sakazakii* as determined using spiral plating and ecometric techniques. *Appl Environ Microbiol.* 2005;71:7661-7669.
12. Teramura H, Fukuda N, Okada Y, Ogihara H. Comparison of chromogenic selective media for the detection of *Cronobacter* spp. (*Enterobacter sakazakii*). *Biocontrol Sci.* 2018;23:27-33.
13. Torres-Chavolla E, Ramírez-Cerda E, Gutiérrez-Rojo R. Isolation and identification of *Enterobacter sakazakii* in infant milk formulas. *Foodborne Pathog Dis.* 2007;4:164-168.
14. Derda R, Lockett MR, Tang SKY, Fuller RC, Maxwell EJ, Breiten B, et al. Filter-based assay for *Escherichia coli* in aqueous samples using bacteriophage-based amplification. *Anal Chem.* 2013;85:7213-7220.
15. Yan W, Qian J, Ge Y, Ye K, Zhou C, Zhang H. Principal component analysis of MALDI-TOF MS of whole-cell foodborne pathogenic bacteria. *Anal Biochem.* 2020;592:113582.